

### DUAL™ Monitor Strip

*A multi-parameter chemical indicator for use with Steam or EO Gas sterilization processes.*

#### Description

The DUAL™ Monitor Strip provides distinct color change when exposed to Steam or EO Gas. Indicator turns dark brown in Steam and orange in EO Gas. The eight inch strip can be used for larger packs, or torn in half at perforation for smaller packs.

#### Features & Benefits

- Multi-parameter ink
- Reduces inventory
- Error free
- Saves money over single process indicators


#### Frequency of Use

"An internal CI should be used within each package, tray or rigid container system to be sterilized. The CI should be placed in that area of the package, tray or containment device considered to be least accessible to steam penetration."

– ANSI/AAMI ST79:2006/A1:2013. *Comprehensive guide to steam sterilization and sterility assurance in health care facilities.*

#### ORDERING INFORMATION

Cat. No.	Description	Packaging
DMS-250	DUAL™ Monitor Strip	250 strips/bx
DMS-250L	DUAL™ Monitor Strip - Laminated	250 strips/bx

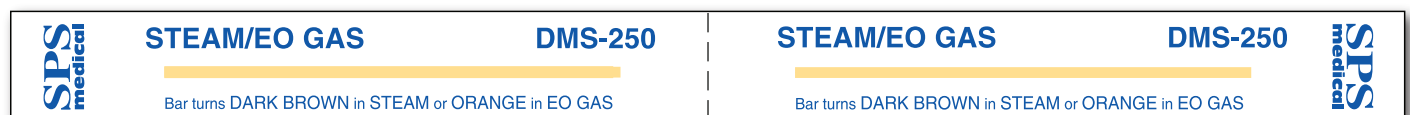

Unprocessed



Processed in **STEAM**, the indicator will turn as dark or darker than this **color standard**.



Processed in **EO GAS**, the indicator will turn as dark or darker than this **color standard**.

DUAL is a trademark of SPSmedical Supply Corp.