

STERILIZATION



Giving You 3x The Confidence!



Sure-Check® Sterilization Pouches



Internal/external multi-variable indicators only change color when all 3 criteria for sterilization have been met:


**Time
Temperature
Steam**

The only way to be SURE about sterilization.

- Built-in multi-parameter internal/external process indicators that change color when all three criteria for sterilization have been met – time, temperature
- Manufactured from medical-grade paper with extra-wide chevron to protect against instrument breaching, thumb-notch for easy opening, 6 corner tack seals to prevent curling and other potential contaminants from entering pouch upon opening and a perforated fold to help ensure the adhesive strip forms an airtight seal
- Eliminates the extra step of adding separate multi-parameter indicators to each pouch
- Once sterilized, contents will remain sterile for up to 3 years, provided the integrity of the pouch has not been compromised
- Meets CDC Guidelines for Disinfection and Sterilization in Healthcare Facilities 2008
- Meets ANSI/AAMI ST79 Comprehensive Guide to Steam Sterilization and Assurance in Healthcare Facilities
- Tested and certified by independent third-party laboratory
- For use with steam and EtO sterilization processes only
- Available in 14 sizes

Sure-Check Sterilization Pouches are not considered hazardous waste and can be disposed of in regular trash (40 CFR 261.24 tested using the TCLP, EPA Publication SW-846).

Patent pending

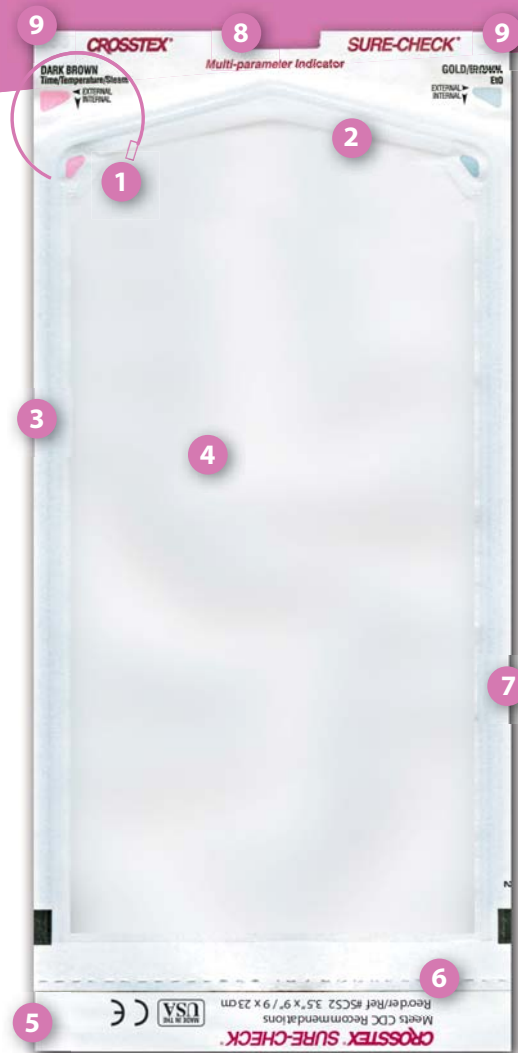

Because You Can't See Sterile®, Be Sure with Sure-Check® Pouches:

- 1 Internal/External Multi-Variable Indicators** - only change color when all 3 criteria for sterilization have been met: **time, temperature and steam.**
- 2 Secure Seal** - Extra wide chevron to protect against instrument breaching.
- 3 Medical Grade Paper** - Wet strength treated, virgin medical grade paper that meets medical pouch standards.
- 4 Transparent Film** - Allows a clear view of the contents as well as the internal indicators.
- 5 Wide Self-Seal Adhesive Strip** - Wider sealing area of strip securely seals pouch and provides sterilization integrity of contents during storage. Contents will remain sterile for three years provided the integrity of the pouch is not compromised.
- 6 Perforated Fold** - Allows for accurate folding of adhesive strip to form an airtight seal, maintaining the integrity of the pouch during storage after sterilization.
- 7 Complete Film/Paper Separation** - The high quality of the materials and manufacturing process assures a consistently clean separation (peeling) of the film from the paper, maintaining the sterile integrity of the contents.
- 8 Thumb Notch** - Assures easy opening of pouch.
- 9 Corner Tack Seals** - Three tack seals on each corner prevent curling and possible separation after sterilization. In addition they prevent the collection of dust and other contaminants in the sealed area.

AVAILABLE IN 14 SIZES!

Description	Quantity	Ref. #
2.25" x 4"/6 x 10 cm	200/Box	SCXX2
2.75" x 9"/7 x 23 cm	200/Box	SCX2
3.5" x 5.25"/9 x 13 cm	200/Box	SCXS2
3.5" x 9"/9 x 23 cm	200/Box	SCS2
3.5" x 9"/9 x 23 cm	500/Box	SCS25
3.5" x 22"/9 x 56 cm	100/Box	SCL3222
5.25" x 6.5"/13 x 17 cm	200/Box	SCW2
5.25" x 10"/13 x 25 cm	200/Box	SCM2
5.25" x 15"/13 x 38 cm	200/Box	SCL5152
7.5" x 13"/19 x 33 cm	200/Box	SCL2
8" x 16"/20 x 41 cm	200/Box	SCL8162
10" x 15"/25 x 38 cm	100/Box	SCL10152
10.5" x 17"/27 x 43 cm	100/Box	SCL10172
12" x 15"/30 x 38 cm	100/Box	SCL12152
12" x 18"/30 x 46 cm	100/Box	SCL12182

Not Made with Natural Rubber Latex

For Questions Relating to Sure-Check® Sterilization Pouches Please Contact Crosstex Customer Care @ 800-819-3336 or Email customer care@crosstex.com

Visit CrosstexBMS.com for Educational Videos

Sure-Check® and Because You Can't See Sterile® are registered trademark of Crosstex International, Inc.